


GETTING STARTED

1. Unpack the contents of the PinkWave™ box:
 (1) PinkWave™ cordless handpiece and battery pack, (1) charging base, (1) charging cord,
 (5) autoclavable light shields, (100) PinkWave™ Barrier Sleeves and (1) User's Manual.
There will also be (3) pairs of protective eyewear included in your shipment.
2. Plug charger into charging base and plug into wall outlet.
 Set PinkWave™ handpiece into the base to charge.
3. The pink ring on the charging base should pulsate, indicating that the light is charging.
 The pink ring will turn solid when PinkWave™ is fully charged.


NOTE: Low battery indicator light will flash red when battery reaches 33% battery life available.

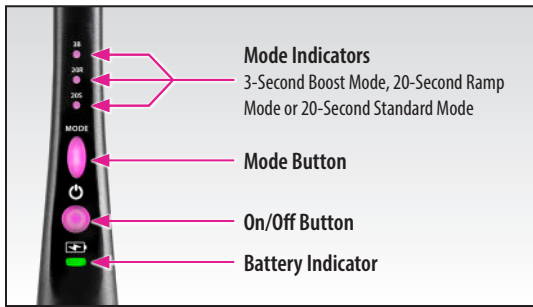

SETTING A MODE

1. Identify the mode button. Use a short press to toggle between settings.
 A light will illuminate the setting that is being selected.

- 3B → 3-second Boost Mode
- 20R → 20-second Ramp Mode
- 20S → 20-second Standard Mode
- All → Transillumintor

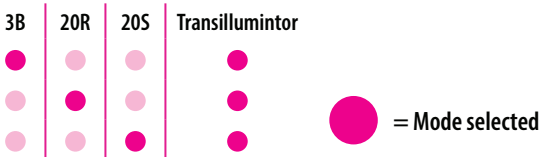

NOTE: Setting is memorized as default until charged.



Refer to User's Manual for detailed instructions.

91047-AP-ENG (5)

A majority of major components are Made in Israel. Final construction performed in USA.

OPERATING LIGHT

1. Once setting is selected, give one short press on either the front or back On/Off button. This will start the curing cycle that was chosen.

The light will automatically turn off when the curing cycle is complete.

On standard and ramp modes, the device will beep once at 5 seconds, twice at 10 seconds and three times at 15 seconds.



2. Press On/Off to power off at anytime.

On/Off Button



NOTE: Make sure light shield is in place before turning on.



CAUTION: DO NOT operate PinkWave™ without the light shield and Vista Apex provided protective eyewear.

TRANSILLUMINATOR

1. To switch to Transilluminator mode, use the mode button to select mode 4, where all lights are illuminated.
2. Press On/Off.

On/Off Button



SAFETY AND DISINFECTION

1. Always have an autoclavable light shield in place before turning device on.
2. Device may be wiped down with a standard disinfecting wipe. Refer to User's Manual for full details and sanitizing protocol.
3. Only PinkWave™ Barrier Sleeves should be used with the PinkWave™ Curing Light.

NOTE: Light shields are autoclavable. Please reference your autoclave instructions for sanitizing protocol.



CAUTION: Never spray disinfectant directly on surface.



CAUTION: DO NOT operate PinkWave™ without the light shield and Vista Apex provided protective eyewear.



Refer to User's Manual for detailed instructions.

91047-AP-ENG (5)

A majority of major components are Made in Israel. Final construction performed in USA.