

Crown Cementation Protocol


VistaApex Solutions Quick Chairside Guide

1. PREP & TRY-IN

1. Prep tooth

2. Place **Stat-Pak™ Knitted Retraction Cord*** if needed3. Use **Quick-Stat™ FREE*** for fluid control

4. Take impression or scan

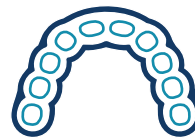

5. Remove temporary > clean prep > try in crown > adjust fit/contacts



*VistaApex Products: **Stat-Pak™ Knitted Retraction Cord** 302611 • **Quick-Stat™ FREE** 504600

2. CONDITION THE RESTORATION BY CROWN TYPE

Zirconia

1. Rinse and dry after try-in
2. Apply **ZR-C™ Universal Cleanser*** for 20 sec
3. Rinse thoroughly and dry
4. Apply **ZR-P™ Universal Primer***
5. Air dry




Porcelain / Etchable Glass Ceramic

1. Confirm material/lab treatment
2. Clean after try-in
3. Apply **Porcelain Etch*** as indicated
4. Rinse and dry
5. Apply silane or **ZR-P™ Universal Primer***
6. Air dry



Metal / PFM

1. Rinse and dry after try-in
2. Apply **ZR-C™ Universal Cleanser*** for 20 sec
3. Rinse thoroughly and dry
4. Apply **ZR-P™ Universal Primer***
5. For PFM, follow exposed intaglio substrate

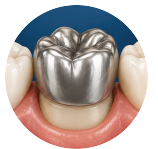

*VistaApex Products: **ZR-C™ Universal Cleanser** 90779 • **ZR-P™ Universal Primer** 90577 • **Porcelain Etch** 502500

3. PREPARE THE TOOTH



Clean / Etch

Use **Velvet Etch™** as indicated.
Clean & Boost® may also be used to clean and condition the prep.



Desensitize

Apply **DesensiMAX™** after etching when indicated.



Optional Adhesive

Apply **RE-GEN™ Universal Adhesive***—doctor preference/IFU dependent.

*VistaApex Products: **Velvet Etch™** 503315 • **Clean & Boost®** 90417 • **DesensiMAX™** 90413 • **RE-GEN™ Universal Adhesive** 90637

4. BIOACTIVE CEMENT

RE-GEN® BIOACTIVE Cement

Set restoration fully, maintain pressure, remove excess, floss contacts as appropriate.

More than a cement, a regenerative solution.

- 90920 Shade A2
- 90921 Shade WO (White Opaque)
- 90922 Shade TR (Translucent)



5. CURE & FINISH



1. Light cure accessible margins as indicated with **PinkWave™ Curing Light***



2. Remove residual cement



3. Recheck contacts and occlusion

*VistaApex Products: **PinkWave™ Curing Light** 90734

