



**Quick-Stat™
Liquid**

Quick-Stat™ Liquid is an aqueous based 15.5% Ferric Sulfate solution, which is excellent for soaking retraction cord, eliminating sulcular fluid contamination, aiding gingival retraction and reducing mucous membrane inflammation. Quick-Stat™ Liquid achieves hemostasis in seconds!
*Compare to Ultradent®: Astringedent®**

- 317120-2 Quick-Stat™ Liquid 30mL Bottle**
• 1 x 30mL Bottle
- 504155-2 Quick-Stat™ Liquid 30mL Bulk Syringe**
• 1 x 30mL Liquid Bulk Syringe



**Stat-Pak™
Knitted
Retraction Cord**

Stat-Pak™ Knitted Retraction Cord is used to retract sulcular tissue, protect tissues and arrest sulcular fluid seepage. With compressible, interlocking chains, knitted Stat-Pak™ Dental Retraction Cord carries more hemostatic agent to the sulcus and remains in-situ, even when touched by a bur or diamond. Stat-Pak™ knitted retraction cord expands after packing, conforming to sulcular anatomy and is memory free for easy placement and precision.
Each bottle contains 100 in (245 cm)

SIZE / COLOR:
000 (Black) 00 (Navy Blue) 0 (Lavender) 1 (Light Blue) 2 (Yellow)

KIT OPTIONS

- 302611-2 #000 Black Stat-Pak™ Cord**
• 100 in (245 cm) 000 - Black Stat-Pak™ Cord
- 302612-2 #00 Navy Blue Stat-Pak™ Cord**
• 100 in (245 cm) 00 - Navy Blue Stat-Pak™ Cord
- 302613-2 #0 Lavender Stat-Pak™ Cord**
• 100 in (245 cm) 00 - Lavender Blue Stat-Pak™ Cord
- 302614-2 #1 Light Blue Stat-Pak™ Cord**
• 100 in (245 cm) 1 - Light Blue Stat-Pak™ Cord
- 302615-2 #2 Yellow Light Blue Stat-Pak™ Cord**
• 100 in (245 cm) 2 - Yellow Stat-Pak™ Cord



**Stat-Pak™
+ E Braided
Retraction Cord**

Stat-Pak™ + E Braided Retraction Cord is 100% cotton braided to prevent strand separation and snagging of packing instruments and is impregnated with stabilized Epinephrine solution. Stat-Pak™ + E is designed to suppress bleeding and control sulcular fluids, while conforming to sulcular anatomy.
Each bottle contains 72 in (176 cm)
Compare to Ultradent®: Ultrapak® E Cord

KIT OPTIONS

- 302633-2 #00 White Stat-Pak™ + E Cord**
• 72 in (176 cm) 00 - White Stat-Pak™ Cord + E
- 302632-2 #0 White Stat-Pak™ + E Cord**
• 72 in (176 cm) 0 - White Stat-Pak™ Cord + E
- 302631-2 #1 White Stat-Pak™ + E Cord**
• 72 in (176 cm) 1 - White Stat-Pak™ Cord + E



**Vista-Pak™
Packers**

Vista-Pak™ Packers have a thin design – ideal for placing knitted retraction cords.
• Plain edge: to place and compress cord
• Serrated edge: for ease of use with knitted cord
Comes in Standard and Circllet Styles
• Circllet style has 3 packing sides for more efficient circular packing of the prep
Compare to Ultradent®: Ultrapak® Packers

- 302703-2 1 x Standard Plain Packer**
- 302704-2 1 x Standard Serrated Packer**
- 302701-2 1 x Circllet Plain Packer**
- 302702-2 1 x Circllet Serrated Packer**

Make Your Patients



...Bleeding




Vista Apex offers a hemostatic solution for every type of bleeder.

| | Quick-Stat™ Gel | Quick-Stat™ Liquid | Quick-Stat™ Free |
|------------------------------|---|---|---|
| Active Ingredients | 20% Ferric Sulfate | 15.5% Ferric Sulfate | 25% Aluminum Magnesium Chloride |
| pH | < 1.0 | < 1.0 | < 1.0 |
| Viscosity | Viscous gel | Aqueous solution | Clear, viscous gel |
| Where to Use | Posterior region | Posterior region | Anterior, aesthetic regions |
| Indications for Use | Gel lubricates cord and allows for precise placement. | Soaking retraction cord | Does not leave behind any residue or interfere with bonding. Ideal for use with self-etch (no rinse) adhesives. |
| Method of Application | Syringe with Stat-Flo™ tip | Syringe with Stat Flo™ tip | Syringe with Stat-Flo™ tip |
| Comparison | Ultradent®ViscoStat™ | Ultradent®Astringent™, PremierHemodent™ | Ultradent®ViscoStat™ Clear,Pascal® Gel Cord® Clear 25%,Dentsply Sirona Hemoban™ Gel |

Technique

STEP 1

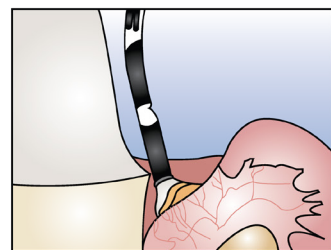


Scrub the hemostatic agent firmly against the sulcus until bleeding is controlled.

STEP 2




A tip with a padded end allows for gentle pressure to be applied, allowing hemostatic agent to penetrate cut capillaries.

STEP 3

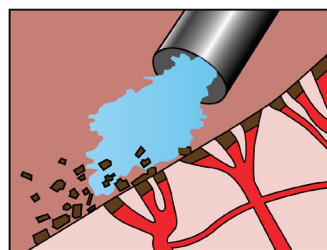

Use air/water to remove coagulum. If bleeding persists, repeat application.

STEP 4

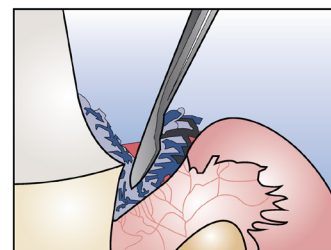

Once complete hemostasis has been achieved, use Stat-Pak™ retraction cord for tissue displacement.



Quick-Stat™ FREE

Quick-Stat™ FREE is the first clear hemostatic agent with shortened hemostasis time. Metal ions in hemostatic agents cause the precipitation of proteins. The positively charged metal ions bind with the negative charge of the blood/serum proteins and you get coagulation because the denatured proteins help accumulate platelets and fibrinogen to speed up clot formation. Not only that, the proprietary surfactant in Quick-Stat™ FREE allows better penetration to more effectively plug cut capillaries. This is the Astringent Hybrid Technology that makes Quick-Stat™ FREE 20% faster than other clear hemostatic agents. For controlling bleeding, especially in the aesthetic zone, Quick-Stat™ FREE is the clear choice. Compare to Ultradent®: ViscoStat® Clear



504600-2

Quick Stat™ FREE Standard Kit

- 4 x 1.2mL Prefilled Syringes
- 8 x Stat-Flo™ Tips



504640-2

Quick-Stat™ FREE Value Pack

- 10 x 1.2mL Prefilled Syringes
- 50 x Stat-Flo™ Tips



504650-2

Quick-Stat™ FREE Bulk Syringe Kit

- 1 x 30mL prefilled syringe,
- 20 x Stat-Flo™ tips
- 20 x 1.2mL empty syringes
- 1 x docking port



504660-2

Quick-Stat™ FREE Bulk Syringe Refill

- 1 x 30mL prefilled syringe
- 2 x docking port



Quick-Stat™ Gel

Our Quick-Stat™ Gel is a highly effective 20% ferric sulfate gel used for immediate and profound control of bleeding. This gel works on both hard and soft tissue and is highly effective for the control of bleeding during impressions, restorations, and crown and bridge procedures. Quick-Stat™ Gel will not flow onto surrounding tissue and quickly controls bleeding. Delivered with Stat-Flo™ applicator tips, this hemostatic agent is easy to use and incredibly effective. Compare to Ultradent®: ViscoStat®



502804-2

Quick-Stat™ Gel Standard Kit

- 4 x 1.2mL Prefilled Syringes
- 8 x Stat-Flo™ tips



502800-2

Quick-Stat™ Gel Value Pack

- 40 x 1.2mL Prefilled Syringes
- 50 x Stat-Flo™ Tips



503550-2

Quick-Stat™ Gel Bulk Syringe Kit

- 1 x 30mL Prefilled Syringe
- 20 x Stat-Flo™ Tips
- 20 x 1.2mL Empty Syringes
- 1 x Docking Port



503555-2

Quick-Stat™ Gel Bulk Syringe Refill Kit

- 1 x 30mL Prefilled Syringe
- 1 x Docking Port