

pink differently




pinkwave™


VISTA | APEX

Patented

Goodbye Blue, Hello Pink

PinkWave was distinctly designed to enhance the efficacy of curing procedures. PinkWave is equipped with Apex's Patented QuadWave™ Technology which leverages four different wavelengths to ensure peak performance. Along with the standard blue light, PinkWave also employs red, near infrared (NIR) and UV light, which together, make up the only pink light on the market. This innovative approach to curing has never been achieved before and results in both doctor and patient benefits.

- Built in Transilluminator
- 3 curing modes — Boost, Standard, and Ramp
- 1640 mW/cm² on Standard mode and 1865 mW/cm² on Boost mode

PinkWave™ Compared to Standard Blue Light



Decrease in
Shrinkage



Increase in
Polymerization



Increase in
Curing Area



Depth of Cure




Lightweight

How you cure matters

quadwave™
Technology

PinkWave™ QuadWave™ Technology leverages four different wavelengths to ensure peak performance:

- UV-photo Initiators
- CQ Initiator
- Enhanced Polymerization
- NIR – Enhanced Polymerization

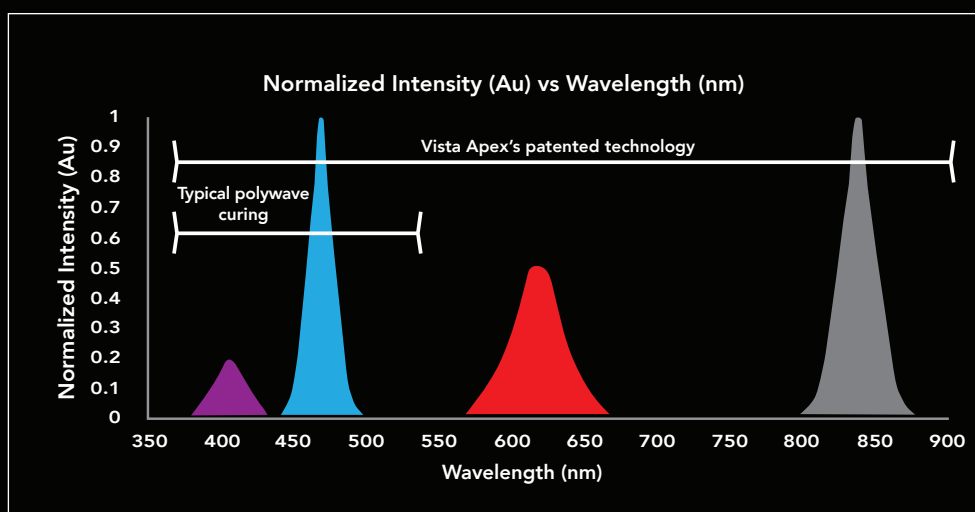

3 Year Warranty

Quadwave™ the future

Decreased Shrinkage


Decreased shrinkage of the materials results in a better restoration with a decreased chance of failure.



Decreased marginal staining, lower chance of microleakage and secondary caries



Reduced patient post-operative sensitivity

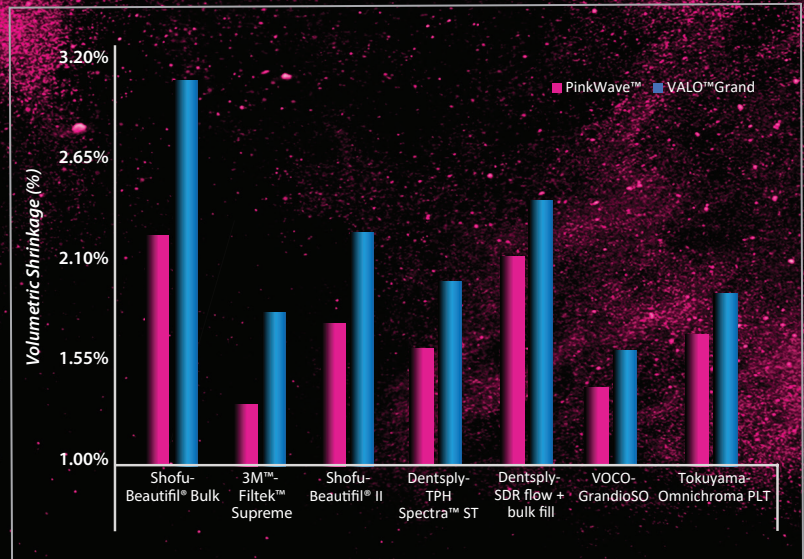

Increased Polymerization


This increase comes from the NIR wavelength emitted and leads to an overall decrease in restoration failure



Harder restoration and greater wear resistance



Less microleakage and water absorption

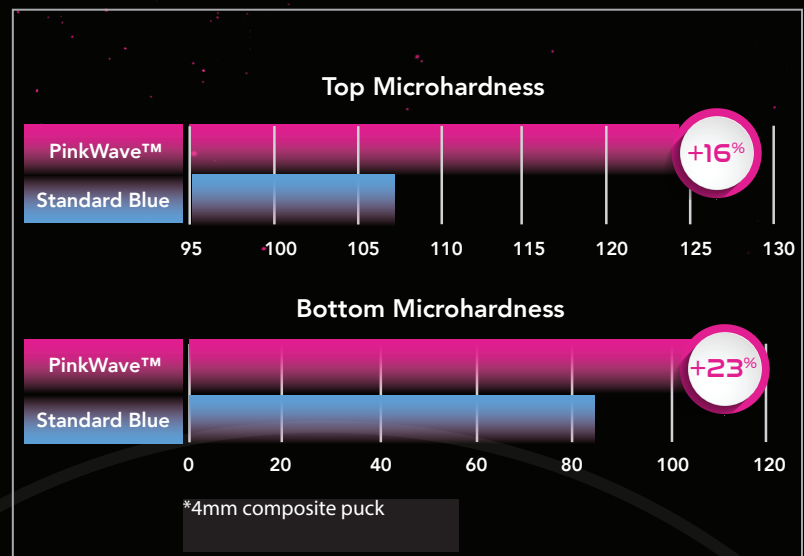

Safe to Use


Blue light output is limited by the ability to severely damage pulpal tissue if too much power is emitted. Since the NIR wavelength in PinkWave is not absorbed by pulpal tissue, PinkWave is able to safely deliver more power.



Patient comfort



Pulp Safety

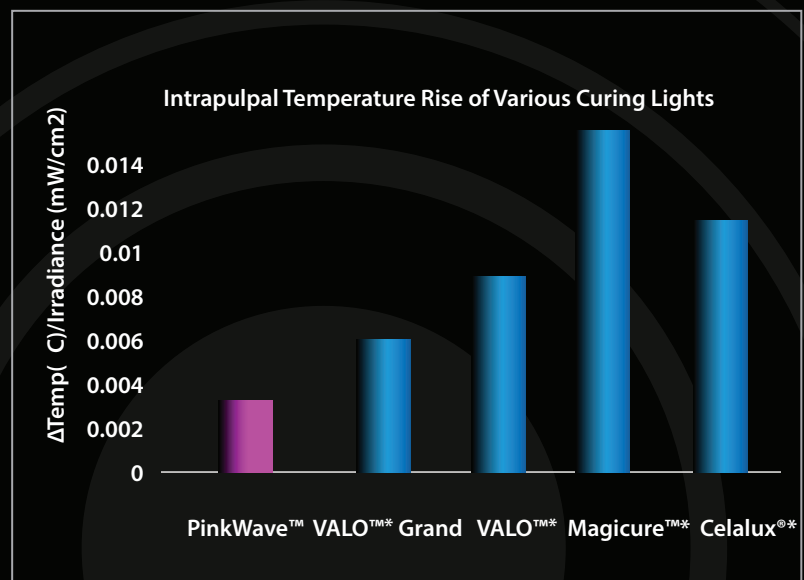

Largest Curing Area



PinkWave™
113 mm²



VALO™* Grand
107 mm²



VALO™*
78 mm²



Valiant™*
44 mm²



Bluephase
41 mm²



SmartLite
38 mm²



Celalux®* 3
36 mm²

Transilluminator



Easy posterior access



Ordering Information

PinkWave™ Kit 90734

- (1) Cordless Curing Light
- (1) Induction Charging Base
- (1) Power Cord
- (5) Autoclavable Light Shields
- (100) Barrier Sleeves



*Percent subject to change on materials used.

*VALO™ is a registered trademark of Ultradent® Products, Inc.. Valiant™ is a trademark of Inter-med, Inc.

*Bluephase® is a registered trademark of Ivoclar Vivadent Inc.

*SmartLite® is a registered trademark of Dentsply Sirona. Celalux® is registered of VOCO GmbH.

*Celalux® is registered of trademark VOCO GmbH.



1-877-APEX 123 (1-877-273-9123)
www.apexdentalmaterials.com