

Prepare for

# IMPACT



IMPACT AIR 45®



## Specifically Engineered for Surgical Procedures

- + Apicoectomy (root canal drainage)
- + Odontoplasty
- + Osteoplasty (bony contouring)
- + Root Resection



VISTA APEX 8774184782 | [vistaapex.com](http://vistaapex.com)

Scan for  
more info >



# IMPACT AIR 45<sup>®</sup> Equipment



## IMPACT AIR 45<sup>®</sup> Highspeed Handpieces

The IMPACT AIR 45<sup>®</sup> has a unique 45° angled head, allowing for maximum access and visibility. It is available with fiberoptics, is fully autoclavable, has an office-replaceable turbine (minimizing downtime) and uses standard friction grip burs (featuring a push-button bur release).

The IMPACT AIR 45<sup>®</sup> is also designed for use in periodontal and endodontic procedures, including Apicoectomy (root canal drainage), Osteoplasty (bony contouring), Odontoplasty and Root Resection.

A water or saline line inside the handpiece handle directs a pure water jet (not a mist) directly onto the cutting area. Air is exhausted through the back of the handpiece.

- Fits A.D.A. Standard Coupling Types
- I.S.O. Type A (3 hole)
- I.S.O. Type B (4 hole or 5 hole fiberoptic)
- Fits our 4 hole, 5 hole, or 6 pin coupler or an equivalent KaVo MULTiflex<sup>®</sup> style coupler
- Compatible with friction grip shank burs (.0627" - .0630" diameter shank)
- Supply Requirements: Dry filtered air or nitrogen
- Operating Pressure: 32psi - 40psi
- Speed: 400,000-500,000 rpm
- Chuck type: Mid-spindle, tooth-jawed grip



### HOSING CONNECTIONS



I.S.O. Type C  
(6 Pin)  
(with adapter)



I.S.O. Type B  
(4 Hole)



I.S.O. Type B  
(5 Hole)

PRODUCT	SKU
#401 IMPACT AIR 45 <sup>®</sup> Highspeed Fiberoptic Handpiece (5 Hole) <i>Includes: (1) IMPACT AIR 45<sup>®</sup> Fiberoptic Handpiece, (1) Head Cap Wrench</i>	DE-401-FO
#403 IMPACT AIR 45 <sup>®</sup> Highspeed Handpiece (4 Hole) <i>Includes: (1) IMPACT AIR 45<sup>®</sup> Highspeed 4 Hole Handpiece, (1) Head Cap Wrench</i>	DE-403-4H
#407 IMPACT AIR 45 <sup>®</sup> Quick Disconnect Highspeed Fiberoptic Handpiece <i>Includes: (1) IMPACT AIR 45<sup>®</sup> Quick Disconnect Highspeed Fiberoptic Handpiece, (1) Head Cap Wrench</i>	DE-407-FO-QD
#408 IMPACT AIR 45 <sup>®</sup> Quick Disconnect Highspeed Non-Optic Handpiece <i>Includes: (1) IMPACT AIR 45<sup>®</sup> Quick Disconnect Highspeed Non-Optic Handpiece, (1) Head Cap Wrench</i>	DE-408-QD

**Vista Apex is the only authorized repair facility for the IMPACT AIR 45<sup>®</sup>.**  
**Contact us at 877.418.4782 for all your handpiece servicing needs.**

**FREE STANDARD SHIPPING ON ORDERS OVER \$400 | CALL TOLL FREE 877.418.4782 OR VISIT [VISTAAPEX.COM](http://VISTAAPEX.COM)**

**ORDER ANYTIME THROUGH YOUR PREFERRED DEALER**

# IMPACT AIR 45® Equipment

## IMPACT AIR 45® Accessories



### Quick Disconnect Couplers

For the IMPACT AIR 45® Quick Disconnect.



PRODUCT	SKU
#627 Quick Disconnect Fiberoptic Coupler (compatible with 6 pin hose)	DE-627COUPLER
#628 Quick Disconnect Fiberoptic Coupler (compatible with 4 hole hose)	DE-628COUPLER
#629 Quick Disconnect Fiberoptic Coupler (compatible with 5 hole hose)	DE-629COUPLER



### Quick Disconnect Parts

For the IMPACT AIR 45® Quick Disconnect.



PRODUCT	SKU
#630 Quick Disconnect Coupler O-Ring	DE-630



### Burs

For all IMPACT AIR 45® handpieces.

\* #520 not shown.



PRODUCT	SKU
#500 Surgical Bur Kit <b>Includes:</b> (1) #4 bur, (1) #6 bur, (1) #701 bur, (1) #702 bur, (1) #1701 bur and (1) #1702 bur	DE-500BUR-KIT
#514 #6 Surgical Length FG (25mm) Bur (5Pk)	DE-514BUR6
#515 #1702 Surgical Length FG (25mm) Bur (5Pk)	DE-515BUR1702
#516 #1703 Surgical Length FG (25mm) Bur (5Pk)	DE-516BUR1703
#517 #702 Surgical Length FG (25mm) Bur (5Pk)	DE-517BUR702
#520 #1701 Surgical Length FG (25mm) Bur (5Pk)	DE-520BUR1701

### Cleaning & Maintenance Kit

For any IMPACT AIR 45® handpieces.



PRODUCT	SKU
#638 IMPACT AIR 45® Cleaning & Maintenance Kit <b>Includes:</b> (1) Spindle Cleaner, (1) Sample Lubricant Bottle, (1) Head Cap Wrench, (1) Blue Fluorosilicone O-Ring, (1) Black Viton O-Ring, (1) 0.005 Teflon Washer, (1) 0.015 Teflon Washer	DE-638


SCAN THE QR CODE TO GO DIRECTLY TO EACH PRODUCT LISTING OR VISIT [VISTAPEX.COM](http://VISTAPEX.COM) TO START YOUR ORDER!

ORDER ANYTIME THROUGH YOUR PREFERRED DEALER



# IMPACT AIR 45® Equipment

## IMPACT AIR 45® Accessories

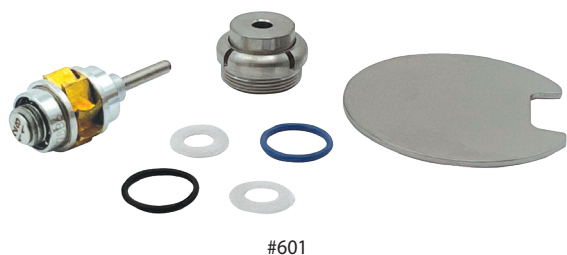

#601




#603



#608

#614

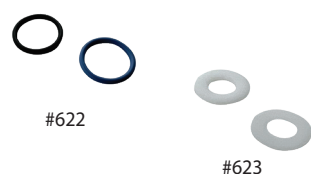

#622

#623




#626



### Parts

For any IMPACT AIR 45® standard handpiece.

PRODUCT	SKU
#601 Model 9 Conversion Kit <b>Includes:</b> (1) M9 Turbine Assembly, (1) Head Cap Assembly, (1) Head Cap Wrench, (1) Blue Fluorosilicone O-Ring, (1) Black Viton O-Ring, (1) 0.005 Teflon Washer, (1) 0.015 Teflon Washer	DE-601CONV
#603 Model 9 Turbine Kit <b>Includes:</b> (1) M9 Turbine Assembly, (1) Blue Fluorosilicone O-Ring, (1) Black Viton O-Ring, (1) 0.005 Teflon Washer, (1) 0.015 Teflon Washer <b>Note:</b> Does not include cap wrench	DE-603TURB
#608 Model 9 Chuck Kit	DE-608CHUCK
#614 Model 9 Head Cap/Wrench Kit <b>Includes:</b> (1) Head Cap Assembly, (1) Head Cap Wrench	DE-614CAP
#622 Turbine O-Ring Kit <b>Includes:</b> (1) Blue Fluorosilicone O-Ring, (1) Black Viton O-Ring	DE-622ORING
#623 Model 9 Teflon Washer Kit <b>Includes:</b> (1) 0.005 Teflon Washer, (1) 0.015 Teflon Washer	DE-623WASHER
#626 Lamp Module For Adapter Kit	DE-626LAMP

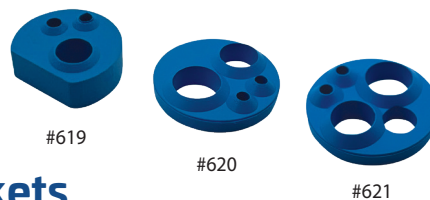

#619

#620

#621

### Gaskets

For any IMPACT AIR 45® standard handpiece.



PRODUCT	SKU
#619 Rear Coupling Gasket 2-Hole Kit	DE-619GASKET2
#620 Rear Coupling Gasket 4-Hole Kit	DE-620GASKET4
#621 Rear Coupling Gasket Fiberoptic Kit	DE-621GASKETFO

### Adapter

For use with the IMPACT AIR 45® or any 4 or 5 hole handpiece on the market.

- Converts your 5 hole style fiberoptic handpieces to the 6-pin fiberoptic hose system
- Compact design
- Easy installation
- Does not affect light quality



PRODUCT	SKU
#625 IMPACT AIR 45® 6-Pin Universal Adapter Kit	DE-625ADAPTER

### Lubricant

For all IMPACT AIR 45® handpieces.

- Penetrates and cleans
- Lubricates
- Protects
- 2,000+ applications per bottle



PRODUCT	SKU
#618 2 Fl. Oz. IMPACT AIR 45® Autoclavable Lubricant	DE-618LUBE



FREE STANDARD SHIPPING ON ORDERS OVER \$400 | CALL TOLL FREE 877.418.4782 OR VISIT [VISTAAPEX.COM](http://VISTAAPEX.COM)

ORDER ANYTIME THROUGH YOUR PREFERRED DEALER