

NEW TO
THE STAGE...




VOYAGER

IRRIGATING TIP

Vista Apex's NEW **Voyager Irrigating Tips** have dual side ports for safety and efficacy. Select gauges are also available in blunt cut. These tips have enhanced flexibility for improved canal navigation, a stainless-steel shoulder for zero flow impedance and maximum visualization, and a Secure-Lock hub. Tips are available in two lengths, 17mm and 27mm, to help fit in limited openings and reach all canal depths.

Our New Rock Star Voyager Tip

- Dual side ports for safety and efficacy, select gauges also available in blunt cut
- Enhanced flexibility to navigate canals
- Available in two lengths to help fit in limited openings reach all canal depths
- Silicone stopper for depth reference
- Pre-bent at the optimal angle for access, easily adjusted by hand or with pliers
- Stainless steel shoulder for zero-flow impedance and maximum visualization
- Secure-Lock hub for reliable syringe attachment

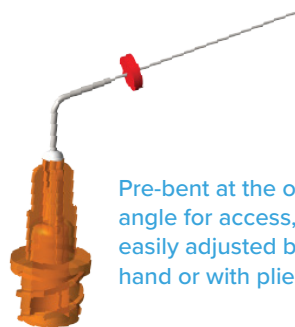

Pre-bent at the optimal angle for access, easily adjusted by hand or with pliers.



SHOULDER BEND



DUAL SIDE PORT

Small flexible cannula to navigate canals and deliver irrigants or sealers.




What makes Voyager Irrigating Tip Better?

| | Available in 4 Gauges | All Gauges Are Dual Side Port | Available in Multiple Lengths | Enhanced Flexibility | Low Cost |
|-----------------------------------|-----------------------|-------------------------------|-------------------------------|----------------------|----------|
| Vista Apex Voyager Irrigating Tip | ✓ | ✓ | ✓ | ✓ | ✓ |
| Ultradent® NaviTip™ | ✗ | ✗ | ✓ | ✓ | ✗ |
| PacDent OptiProbe™ | ✗ | ✗ | ✓ | ✗ | ✓ |
| Benco Dental™ Irrigator Probe | ✗ | ✗ | ✓ | ✗ | ✓ |

Ultradent® and NaviTip™ are registered trademarks of ULTRADENT® Products, Inc. OptiProbe™ is a registered trademark of PacDent, Inc. Benco Dental™ is a trademark of Benco Dental Supply Company.

Spotlight Performance

To ensure safe delivery, this tip was engineered so that irrigants are expressed out the sides of the tip. Using a metal shoulder in place of a large plastic hub, the two lengths of Voyager Tips offer the same full range of working lengths compared to competitors without having to stock as many variations in your operatory. The reduced bulk of the hub also allows for easier access in smaller oral openings.

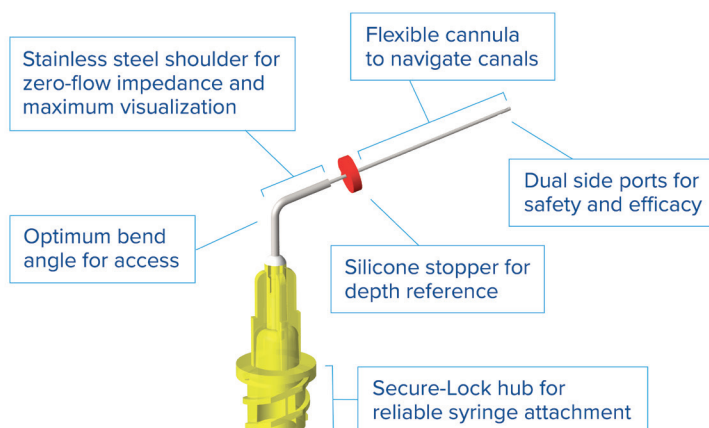

| GA | STYLE | LENGTH | COLOR | PART # |
|------|----------------|--------------|--------|----------------------|
| 27GA | Dual Side Port | 17mm 27mm | Yellow | 315727-2 315827-2 |
| 29GA | Blunt End | 17mm 27mm | Orange | 315929-2 316929-2 |
| 29GA | Dual Side Port | 17mm 27mm | Orange | 315729-2 315829-2 |
| 30GA | Dual Side Port | 17mm 27mm | Purple | 315730-2 315830-2 |
| 31GA | Dual Side Port | 17mm 27mm | Green | 315731-2 315831-2 |



**SCAN
HERE**
to Learn More

TO ORDER:

Call toll free 1-877-418-4782,
visit vistaapex.com/voyager,
or scan the QR code

FAQ

Which tips should I use to deliver NaOCl or EDTA?

Any of the dual side port tips will work depending on your prep size and gauge preference.

Which tips should I use to deliver paste or sealer materials?

27ga and 29ga in dual side port or blunt cut, will be best.

For an even better delivery tip for these materials, check out the MST Visco Tips.

What do the tip colors mean?

Vista Apex tips are color coded by needle gauge for easy identification.

Can Voyager tips be used for Evacuation?

Absolutely, any tips will work but higher suction flow rates will be achieved with the larger 29 and 27 gauges. For even higher flow rates, check out Vista Apex's Micro-flow, Flexi-glide, and Elasti-vac tips.