

Light Curing, Redefined.

PinkWave™



PinkWave was specifically designed to enhance the efficacy of curing procedures. Our patented QuadWave™ technology leverages four different wavelengths (UV, blue, red and, near-infrared) to ensure peak performance. This innovative approach to curing results in benefits for you and your patient.

- 3 curing modes:
Boost, Standard and Ramp
- 1515 mW/cm² on standard mode
and 1720 mW/cm² on boost mode
- Built-in Transilluminator

Depth
of Cure
Up to 8.5mm



VistaApex
SOLUTIONS
Value. Added.


QuadWave™ – The Differentiator

QuadWave™ provides more efficient polymerization and less composite shrinkage. By incorporating near-Infrared light, PinkWave™ shows superior performance across a variety of composite materials.^{1,2} Most curing lights deliver maximum power at the surface of the curing light lens which can be clinically challenging for deeper restorations. With this in mind, PinkWave™ was designed to peak at 4-5mm away from its lens to deliver the power you need deeper into the restoration, regardless of how consistent the distance from the restoration. This results in a harder restoration with greater wear resistance, less microleakage and water absorption.^{1,2}

QuadWave Technology



- UV-photoinitiators
- CQ photoinitiator
- Enhanced polymerization¹
- Red and near-Infrared to minimize shrinkage and increase depth of cure^{1,2}

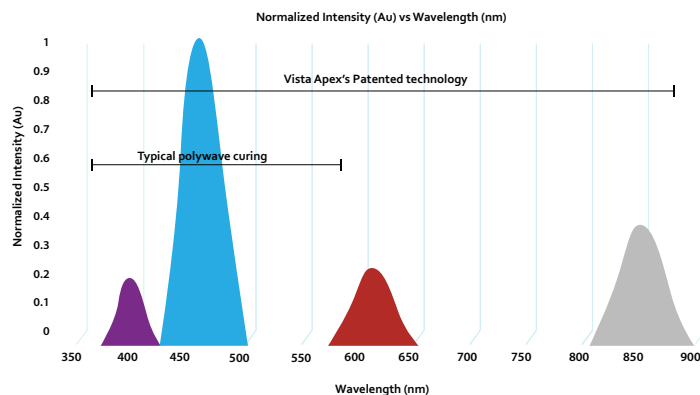

Large
Curing
Area
12.1mm

PinkWave Key Features:

- Depth of cure - up to 8.5mm
- Transillumination - helps detect interproximal cracks
- Safe and lightweight - weighs 88gm and has very low heat emission
- Decreased shrinkage - better restoration and less chance of failure¹
- Low profile head - for easy access to posterior restorations
- Reduced patient post-operative sensitivity



Ordering Information

PinkWave Kit 90734

- (1) Cordless curing light
- (1) Induction charging base
- (1) Power cord
- (5) Autoclavable light shields
- (100) Barrier sleeves
- (3) Pair of eyewear
- (5) Anti-Roll Brackets

VistaApex Solutions

2200 South Street
Racine, WI 53404

Phone: (877) 857-8213

Email: info@vista-dental.com

Website: <https://vistaapex.com>

References

1. Data on File. VistaApex Solutions, Racine WI
2. Ghaffari S, et al. (2025). Dent Mater. 41(5):503-512.

Scan to
learn more

